

XENON²


QUICKIE[®]
Live without limits.



**SUPERIOR
FOLDING
AND
DRIVING
PERFORMANCE**

XENON² – THE STYLISH, HIGH-PERFORMANCE SERIES

Whether you're looking for the very best in driving dynamics or the most compact folding dimensions, there's a Xenon² just for you.

The New Xenon² is available in three different versions:

- The streamlined Xenon² fixed front
- The reinforced Xenon² hybrid fixed front
- The portable Xenon² swing away

Whichever one of the three feature packed models you choose, the sophisticated elegance of the Xenon² will always make you the centre of attention.





XENON² FF THE LIGHTEST

The Lightest folding wheelchair in the world. Open fixed front frame with a cool, clean and streamlined design.



XENON² H (HYBRID) THE STRONGEST

The hybrid is reinforced and has a fixed front frame, no rigidity is compromised and it also has a maximum user weight of 140 kg.



XENON² SA THE ALLROUNDER

The Xenon² SA has the most compact folded dimensions for easy portability by car. The frame is reinforced and has a swing away front, for easier standing transfers in and out of your wheelchair.





THE LIGHTEST FOLDING CHAIR IN THE WORLD...

How do you achieve a folding wheelchair that weighs as little as 8.8 kg? By using premium frame materials like aerospace-grade aluminium and an innovative process that's unique only to Quickie – ShapeLoc Technology.

By heat-treating the aluminium, we're able to increase its strength, which allows for the frame's tube walls to be thinner. Combine this with an oval design and hydroformed ShapeLoc Technology this gives us a lighter and stronger wheelchair than anyone else. Just feel it...

...WITH UNPARALLELED FOLDING PERFORMANCE

At the core of the Xenon² is a new and unique cross-brace that fits so neatly under the seat it's barely noticeable. This creates a very minimalistic, open-frame feel that you'd usually associate with a rigid wheelchair.

It doesn't just end there! The new axle stem has been completely re-designed and combined with the new cross-brace, provides the stiff driving performance and feel of a rigid frame.

Easy to use folding mechanism

The new aluminium latches on the cross-brace will also keep your wheels in parallel when folding, so your Xenon² remains stable no matter what the situation.






STYLE IN MOTION

Personalise your Xenon² and stand out from the crowd with a versatile palette of 32 frame colours. There's also 5 distinctive anodised styles to choose from for the castor fork and castorwheel. You can even accentuate your colour scheme further by matching or contrasting the upholstery trim colour of the EXO backrest upholstery.

→ Personalise your Xenon²
With a selection of 32
eye-catching colours

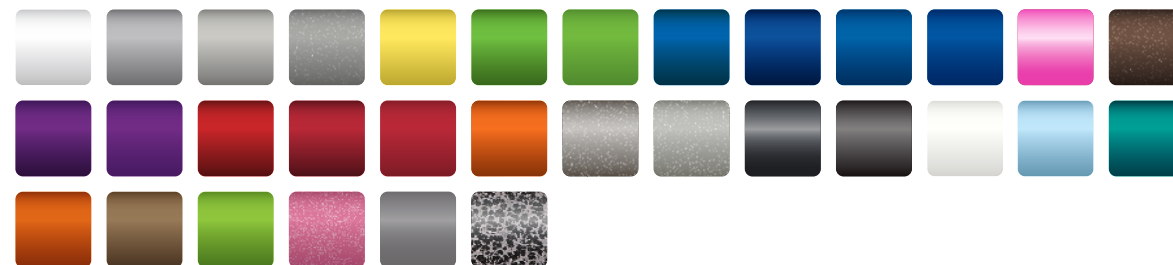



→ **Add some colour to your back upholstery**
This very comfortable and stylish Exo back upholstery can be chosen with 5 different colour trims.



→ **Choose your own castor and fork style**
The anodised front castors and the forks are available in silver, blue, red, orange or black for a unique look.



ADJUSTABLE TO SUIT YOU

Exactly as it should be – designed to be fine-tuned and adaptable to your changing needs.

Do you want to change your chair driving characteristic from passive to more active or the opposite way? Just change the centre of gravity setting by simple repositioning of the axle stem: and don't forget to adjust the back angle for more ergonomic posture.

- ❶ **Centre of Gravity**
Enables a driving characteristic from very active to passive.
- ❷ **Backrest Angle**
The right back angle adjustment is important for your ergonomic posture.
- ❸ **Rear Seat Height**
Exactly as it should be - adjust the height to your personal need.





XENON² FF
THE LIGHTEST



XENON² H (HYBRID)
THE STRONGEST



XENON² SA
THE ALLROUNDER

XENON²

Ultra lightweight wheelchairs
with outstanding folding and
driving performance



Options and features of the Xenon² series – Configure your wheelchair to your individual needs. For further information on product options please check the Xenon² and Xenon² SA order form.



// FOLD-DOWN BACKREST

This half-folding backrest folds to the front and achieves a small package size, making it easier to transport.



// SURGE® LT HANDRIMS

Its Gripton rubber strip on top of the ergonomic oval shaped rim offers maximum traction.



// ANGLE ADJUSTABLE BACK

Stylish designer backrest solution integrated into the back of the frame. Backrest angle from 75° to 103° in 4° increments.



// DESK SIDEGUARD

This flip-up and removable sideguard can be ordered with fixed or height adjustable arm pads.



// PADDED TUBULAR ARMREST

These armrests give you the support when needed and do not hinder during propulsion, they can also swing away.



// CARBON SIDEGUARD

Very strong and lightweight sideguard with a unique look.



// AUTO-FOLDING FOODBORD

This foodplate is angle and depth adjustable and available with carbon or composite. Only available for fixed front frame versions.



// TIE-DOWN BRACKETS FOR TRANSPORTATION

With this option the wheelchair is approved for transportation in a vehicle (crash tested according to ISO 7176-19).

TECHNICAL DATA

	Xenon ² FF	Xenon ² H (Hybrid)	Xenon ² SA
Max. user weight:	110 kg	140 kg	125 kg
Seat width:	32 – 46 cm	32 – 50 cm	32 – 50 cm
Seat depth:	34 - 50 cm	34 - 50 cm	34 - 50 cm
Seat height front:	43 – 55 cm	43 – 55 cm	43 – 56 cm
Seat height rear:	37 – 50 cm	37 – 50 cm	37 – 50 cm
Backrest height:	25 – 47.5 cm	25 – 47.5 cm	25 – 47.5 cm
Backrest angle:	-15° to +13° increments	In 4° increments	In 4° increments
Centre of gravity (COG):	4 – 13 cm	4 – 13 cm	4 – 13 cm
Camber:	0° / 2° / 4°	0° / 2° / 4°	0° / 2° / 4°
Frame Version:	Open frame with fixed front	Reinforced frame with fixed front	Reinforced frame with swing-away front
Chair weight:	From 8.8 kg	From 9.5 kg	From 10.3 kg
Frame colours:	32	32	32
Crash tested:	Yes	Yes	Yes

For further information on the full specification, options and accessories, please check the order form on our website. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have. If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.co.uk or alternatively it is available on request in large text.



Sunrise Medical Limited

Thorns Road, Brierley Hill, West Midlands, DY5 2LD
 Tel: +44 (0) 845 605 66 88 Fax: +44(0) 845 605 55 89
 E-Mail: enquiries@sunmed.co.uk
www.SunriseMedical.co.uk