



Freedom of expression

The *Invacare Spectra XTR2* has been developed to combine a powerful driving performance with compactness and transportability. The remarkable chair has a sophisticated spring suspension designed to provide a smoother, easier ride when outdoors. Its compact design also makes it ideal for indoor use while the easy to remove battery system makes the *Spectra XTR2* extremely quick and simple to service. Key features include:

- Excellent traction
- Compact design
- Highly adjustable seating
- Simply SMART LiNX control system



LiNX®



Yes, you can.®



Sophisticated spring suspension system

The **Spectra XTR2** has a sophisticated spring suspension system to provide outstanding outdoor capabilities and a comfortable drive. It is available with two motor versions for driving speeds of up to 6 or 10 km/h. The remarkable chair also offers a 300mm seat lifter option including powered tilt function up to 30°.

Compact design

Excellent indoor mobility is obtained due to the particularly compact dimensions of the **Spectra XTR2**. When choosing smaller castors, this feature is optimised, providing a width less than 590mm.

Removable seat unit and batteries

With a clever and robust mechanism for detaching the seat from the wheelbase, transportation, storage and services are rendered refreshingly trouble-free. To make servicing and transportation even easier, the batteries with an automatic connecting plug, slide out smoothly and without effort.

Simply SMART LiNX control system

Invacare **LiNX** is our new insight inspired control system that learns over time to bring a revolutionary driving experience. Thanks to its unique modular design, this futureproof system can be seamlessly expanded and tailored to complement changing needs.

The Simply SMART control system allows programming and maintenance updates to be made wirelessly in real time via bluetooth - making the **Spectra XTR2** one of the easiest powerchairs on the market to service and maintain. With easy access to statistical and diagnostic information, **LiNX** allows complete visibility for professionals. Features include:



Icon based menu

The design of the new interface focuses on ease of use, simplicity and reliability. It is very straightforward for everyone to use, including those with limited cognitive abilities.



Total control

For more complex control needs, **LiNX** improves functionality, driving control and performance, while **G-Trac*** (optional) enhances this and ensures the directions indicated by the control are followed by the chair.

FEATURES & OPTIONS



Seat lifter

A seat lifter of 300mm including powered tilt function of up to 30° increases the use of the **Spectra XTR2**. Featuring lower seat to floor heights in conjunction with our **Modulite™** seat system.



Basic seating

With integrated arm and legrest receivers and tension adjustable back support.



Telescopic seat frame

As the telescopic seat frame and the Flex3 backrest plates adjust in different width from 380 to 530mm.



Removable seat unit

After pulling the strap, tilt the seat and disconnect the cables in order to take the seat off. It couldn't be easier!



Removable batteries

Separate battery boxes slide in and out without any effort. High-end sockets connect automatically.



Luggage rack

The very handy luggage rack helps to carry your belongings. It can be mounted in different positions and even removed when not needed.



Swing out anti tippers

A useful feature for storing the wheelbase either in the home or warehouse.



Spring for suspension

A comfortable drive and safe traction are ensured even at high speeds.



Different applications

Choice of two castor sizes depending on customer requirements.

Technical data

Basic seating 440 - 480 mm ¹	410 - 510 mm	450 mm ³	480 mm	275 - 340 mm	290 - 470 mm ^{3,4}	0 - 20°
Modulite™ telescopic seat	380 - 530 mm ²	410 - 510 mm	450 mm ³ / 450 - 750 mm ^{3,10}	480, 540, 560 mm ^{3,11}	245 - 310 mm	290 - 470 mm ^{3,4} 0 - 20° / 0 - 30° ¹⁰
Basic seating	90° - 120°	590 / 620 mm ⁵	930 mm	1160 mm ⁶ 1210 mm ¹⁰	820 / 840 mm ⁵	98 kg ⁹
Modulite™ telescopic seat	90° - 120°	590 / 620 mm ⁵	970 mm	1160 mm ⁶ 1210 mm ¹⁰	820 / 840 mm ⁵	94 kg ⁹
ca. 50 Ah	1600 mm 1630 mm ¹⁰	60 mm ⁷	18% (10°)	2 x 230 W	up to 30 km ⁸	6 / 10 km/h


- 1/ Seat cushion 430 mm
- 2/ Additional 50 mm by armrest adjustment
- 3/ Measured without seat cushion / upholstery
- 4/ Certain legrest lengths are restricted to a special legrest type
- 5/ Width of wheelbase depending on castor size
- 6/ Incl. standard legrests
- 7/ 100 mm with kerb climber
- 8/ According to ISO 7176-4:2008
- 9/ Total weight is depending on configuration
- 10/ With seat lifter or 30° COG tilt module
- 11/ Depending on back support choice
- 12/ **Invacare Modulite™** seat components are cross compatible across the Invacare power wheelchair range and partly with the Invacare Rea[®] family.
- 13/ Compliant with the directive on the Restriction of the use of certain Hazardous Substances in electrical and electronic equipment - RoHS (2011/65/EC).
- 14/ Compliant with the regulation concerning the Registration, Evaluation, Authorization and Restriction of Chemicals.



matrix
SEATING SERIES

Invacare® **Modulite™**

Frame colours

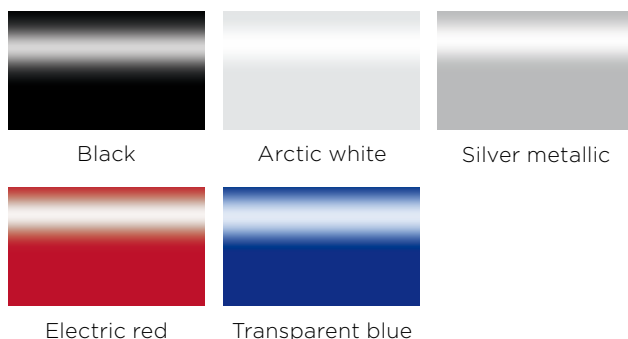

The LiNX word, mark and logos are registered trademarks owned by Dynamic Controls and any use of such marks by Invacare is under license.

For more information about this wheelchair, please see www.invacare.co.uk. For information about typical description, intended use, adjustments and transportation, please refer to the user guide.

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ -
United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. The images in this brochure are for illustration purposes only and may vary slightly from the exact product available.

Spectra XTR2 - UK - 06/2016 - 1568223



Follow us on Twitter:
[@Think_Mobility](https://twitter.com/Think_Mobility)

Follow us on YouTube:
Invacare Ltd - UK



Like us on Facebook:
Invacare UK



Yes, you can.®