



QUICKIE
Live without limits.

ATTITUDE
PURE FREEDOM



EXPERIENCE FREEDOM LIKE NEVER BEFORE

- **Life without limits with the right Attitude** - choose a manual, hybrid or power Attitude add-on hand bike.
- **Roll away with the right Attitude** - lightweight hardware.
- **Connect with the right Attitude** - simple and secure docking.
- **Attitude with style** - 31 different colours.



WHAT'S YOUR ATTITUDE?

Whether you want to exercise with the manual add-on, enjoy the benefits of rehab and intermittent exercise with the power assist or cruise anywhere with the full power add-on, get the right Attitude from Quickie - opening up so many possibilities to live life to the limits. When you ride with Attitude you'll wonder why you've never tried it before!

ATTITUDE Power

FREEDOM TO RIDE

Providing you with 100% powered propulsion so you can hit off-road terrain or cruise the streets in speed for pure enjoyment. Live life to the limits and go anywhere with the Attitude power.



ATTITUDE Manual

FREEDOM TO EXERCISE

Enjoy the benefits of improving your strength, stamina and cardio through recreational cycling with the Attitude manual add-on handbike. Lightweight and easily manoeuvrable, it's the perfect complement to your active wheelchair.



ATTITUDE Hybrid

FREEDOM TO ENJOY

Combine all of the health benefits of manual cycling with the reliability of powered propulsion. Helping to reduce fatigue, the powered motor can be engaged whenever you encounter a tricky hill-climb or just to provide you with that little extra speed when you need it.





ROLL AWAY WITH THE RIGHT ATTITUDE

Speed - Simplicity - Unrivalled Rigidity.

With Quickie compromise is not an option - so we say no to heavy mounting hardware that reduces your everyday performance. Un-dock in seconds and roll away with only two lightweight pins attached to your wheelchair, so whether you're riding with an Attitude handbike or your everyday Quickie wheelchair you're experiencing the very best lightweight performance.




After de-docking only two small pins (170g) remain on your wheelchair.

CONNECT WITH THE RIGHT ATTITUDE

The docking mechanism on the Attitude bike is so simple, it just locks securely to the front of your wheelchair – without any need for parts or cluttered wiring around your knees. So there's nothing coming between you and the comfort that you've come to expect from your wheelchair.

- 1 Lock in the docking pins.
- 2 Swing the upper docking connector downwards.
- 3 Lock the connection by pulling the red lever.

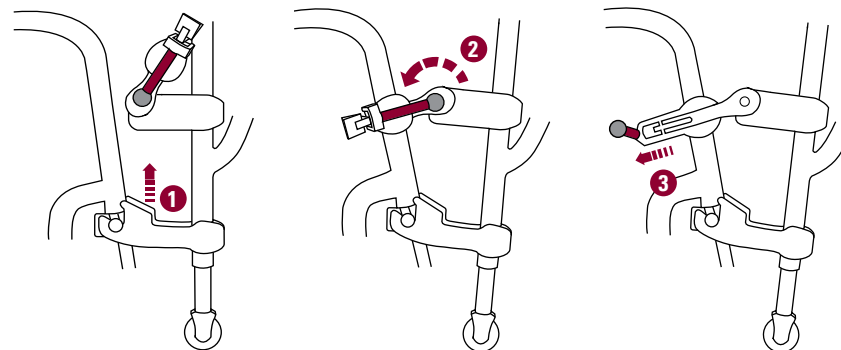

Check out the Quickie Attitude docking video at:

www.sunrisemedical.co.uk/attitudevideo



→ Ready to run the Attitude docking video? Check it out!



MORE FREEDOM

The freedom you'll enjoy with the Attitude bike not only applies to when you're biking. The ergonomic handle and support wheels will help you to manoeuvre the Attitude bike easily with just one hand. So whether you're de-docking your bike or storing it away, it's effortless no matter the situation.



ATTITUDE WITH STYLE

There's a change in attitude with add-on hand bikes! Ultra lightweight to match the performance you'd expect from a Quickie wheelchair, sleek design and 31 head turning colours to match your existing wheelchair or create a look that's truly individual. Hit the streets in style with Quickie Attitude.



ATTITUDE

MANUAL – FREEDOM TO
EXERCISE



Whether you want to get fit or feel more energised - the Manual Attitude add-on gives you the opportunity to improve your health and push your stamina to the limits. Achieve your fitness goals and experience the exhilaration with the Manual Attitude add-on hand bike.

Manual



// DUAL DRIVE (3x10)

30 speed gearing is perfect for all terrains and the disc brakes allow for controlled braking.



// GRIP SHIFT

The easy one hand operation grip shift allows quick and easy gear shifting.



// BOTTLE HOLDER

This bottle holder can hold all standard bottles (available for all Attitude add-on hand bikes).



// PARKING BRAKE

Easy operating parking brake with lock in function as standard on all Attitude add-on hand bikes.

ATTITUDE

HYBRID – FREEDOM TO
ENJOY



Break through your boundaries and explore more than ever before using the Hybrid Attitude. The automatic motor smoothly engages when you start to cycle and enables you to keep up on faster bike rides, tackle rough terrain and steep inclines not possible with the manual version.

Hybrid



// **LITHIUM-ION BATTERY**

14.5Ah Lithium-Ion battery for maximum range performance. The display also shows the battery capacity remaining.



// **ERGONOMIC V-CRANK**

The ergonomic handle position prevents contact with your body particularly when manoeuvring around corners.



// **STRONG 250 WATT MOTOR**

The powerful 250 watt motor supports up to a 23kph top speed, ideal for hilly terrain or longer distances.




// **USEFUL LOW RIDER CARRIER**

The low rider carrier can be used for any available standard saddle bags. This option is available for all Attitude add-on hand bikes.

ATTITUDE

POWER – FREEDOM TO
RIDE

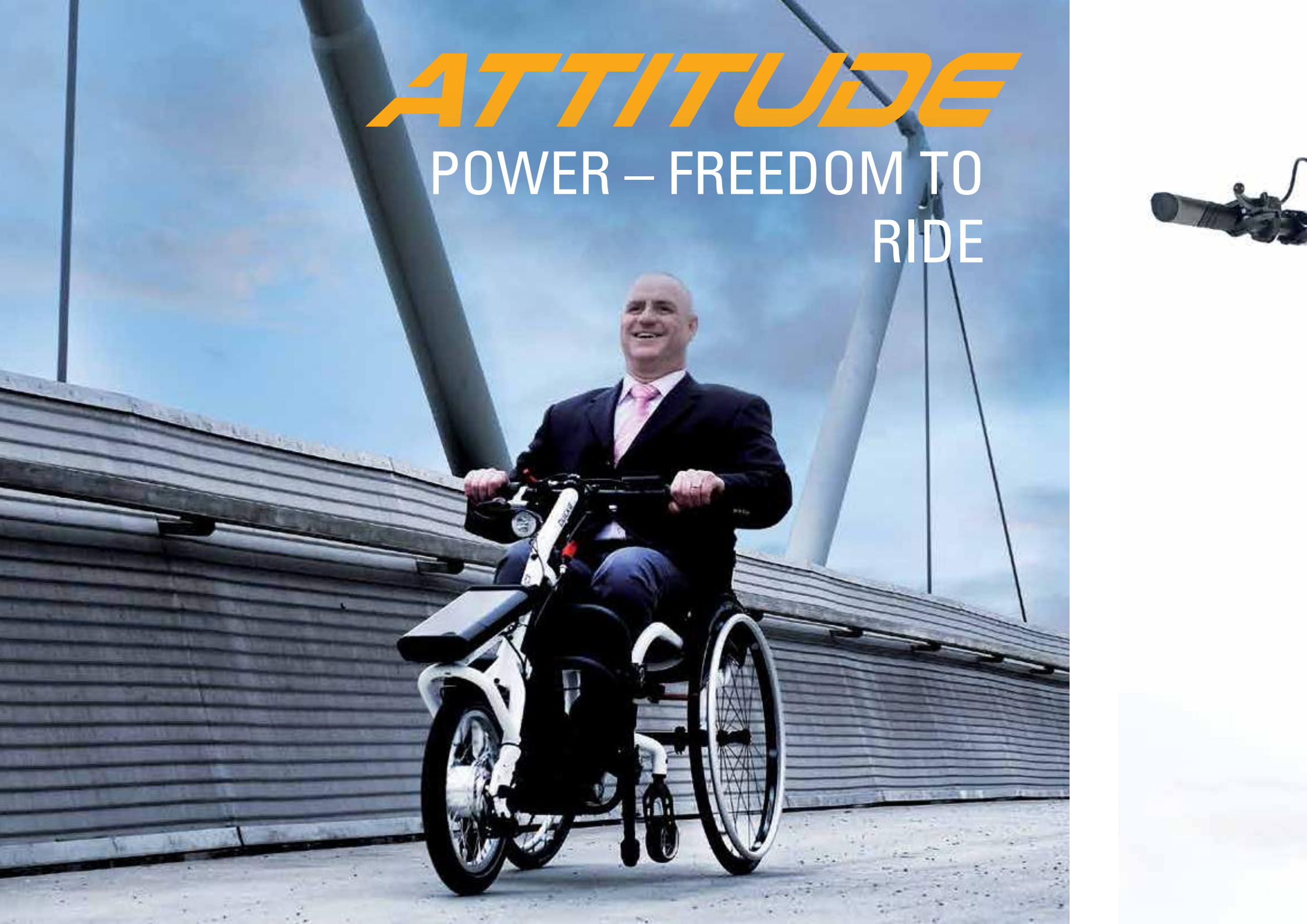

Hit off-road terrain or experience the open road at exhilarating speeds of up to 23kph. Commutes to work, cruising the shops and days out with the family become effortless with the full power add-on.

Power



// STRONG 250 WATT MOTOR

Powerful 250 watt motor with 23kph top speed.



// ADJUSTABLE HANDLE-BAR

The ergonomic handle bar is height adjustable via the stem. The display shows a range of driving options.



// LITHIUM-ION BATTERY

9Ah or 13.4Ah Lithium-Ion battery for large driving range.



// INTEGRATED FRONT LIGHT

The bright front light is powered directly from the battery.

TECHNICAL DATA 46

<i>ATTITUDE</i>	<i>Manual</i>	<i>Hybrid</i>	<i>Power</i>
Max. user weight	100kg	100kg	100kg
Front frame width (adj. + 2cm)	26 – 46cm	26 – 46cm	26 – 46cm
Hardware weight on chair	0.2kg	0.2kg	0.2kg
Total bike weight	14.5kg	21.5kg	15.3kg
Wheel size	20"	20"	16"
Max. motor speed	–	23kph	6kph & 23kph
Motor	–	250W	250W
Battery	–	14.5Ah	9Ah & 13.4Ah
Gearing shifters	Trigger / Grip shift	Trigger / Grip shift	–
Gearing speeds	8 & 30	10	–
V-crank width	470mm	470mm	–
Crank length	170mm	170mm	–
Handle bar width	–	–	570mm
Brake	Disc (160mm)	Disc (160mm)	Disc (160mm)
Parking brake	Yes	Yes	Yes
Frame colour	31	31	31

For further information on product options, please check the order form on our website.
 You can view this document also in a larger font as a PDF file at www.SunriseMedical.co.uk



www.facebook.com/QuickieWheelchairs



www.youtube.com/quickiewheelchairseu



Sunrise Medical Limited
 Thorns Road, Brierley Hill, West Midlands, DY5 2LD
 Tel.: +44 (0) 845 605 66 88 Fax: +44 (0) 845 605 66 89
 E-mail: enquiries@sunmed.co.uk
www.SunriseMedical.co.uk